

Product Data Sheet

Anti-PIEZO1 Antibody

Catalog # Source Reactivity Applications

CQA2997 Rabbit H WB, IF/IC

Description Rabbit polyclonal antibody to PIEZO1

Immunogen Recombinant full length protein of human PIEZO1

Purification The antibody was purified by immunogen affinity chromatography.

Specificity Recognizes endogenous levels of PIEZO1 protein.

Clonality Polyclonal

Conjugation

Form Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol,

and 0.01% sodium azide.

Dilution WB (1/500 - 1/1000), IF/IC (1/50 - 1/200)

Gene Symbol PIEZO1

Alternative Names FAM38A; KIAA0233; Piezo-type mechanosensitive ion channel component 1;

Membrane protein induced by beta-amyloid treatment; Mib; Protein FAM38A

Entrez Gene 9780 (Human)

SwissProt Q92508 (Human)

Storage/Stability Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid


freeze/thaw cycles.

Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC- Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish


COHESION BIOSCIENCES LIMITED



Product Data Sheet

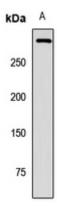

Western blot analysis of PIEZO1 expression in MCF7 (A) whole cell lysates. (Predicted band size: 286 kD; Observed band size: 287 kD)

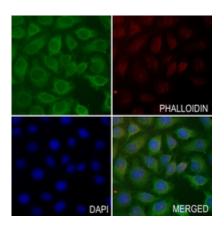

Immunofluorescent analysis of PIEZO1 staining in A549 cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a hidified chamber. Cells were washed with PBST and incubated with a AF488-conjugated secondary antibody (green) in PBS at room temperature in the dark. Phalloidin - AF594 was used to stain Actin filaments (red). DAPI was used to stain the cell nuclei (blue).

Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC-Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference

Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish